


Sterile and adhesive surgical incise drape

# Euroderm incise iodine

Euroderm Incise Iodine respects the same indications and has the same properties of Euroderm Incise. Euroderm Incise Iodine provides continuous antimicrobial broad-spectrum action by the gradual release of iodine molecules for over 7 hours of application. The film is liquid and bacteria proof (also in dry conditions) but it is permeable to gases and water vapour, protecting the area from contamination and leaves the skin to breathe, thereby preventing the accumulation of moisture under the drape. A critical aspect of all surgical interventions is the prevention of infections; Euroderm Incise Iodine reduces the entry of microorganisms into the wound and prevents lateral migration of bacteria into the wound itself. Studies on the effectiveness of Euroderm Incise Iodine on *Staphylococcus Aureus* (11000 CFU/ml) of *Pseudomonas Aeruginosa* (10200 CFU/ml) and *Candida Albicans* (11000 CFU/ml), showed a total reduction of bacteria. Euroderm Incise Iodine is latex free.

## Features and benefits

- Adheres to the edge of the wound.
- Avoid lateral migration of bacteria into the wound.
- Non-reflective surface that provides perfect visibility of the incision area.
- Ideal in cases where the manipulation is required during surgery (hip prosthesis).
- Liquids and bacteria proof.
- Wide range of sizes.

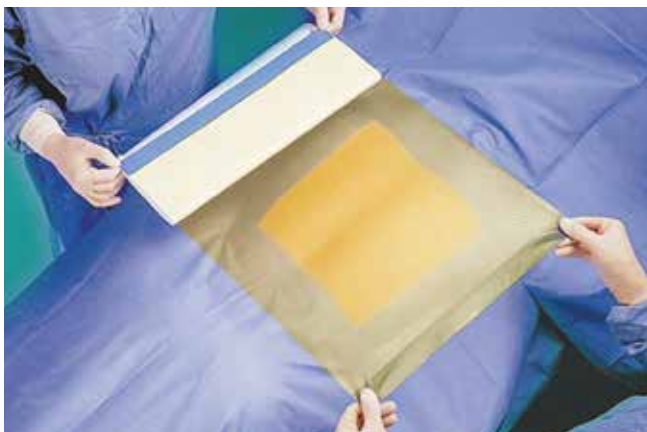

item no.	size	→ adhesive area	pkg.
717910	cm 10 x cm 24	cm 10 x cm 18	10 pz/pcs
717915	cm 15 x cm 26	cm 15 x cm 20	10 pz/pcs
717930	cm 30 x cm 26	cm 30 x cm 20	10 pz/pcs
717940	cm 40 x cm 34	cm 40 x cm 28	10 pz/pcs
717942	cm 42 x cm 42	cm 42 x cm 36	10 pz/pcs
717944	cm 44 x cm 26	cm 44 x cm 20	10 pz/pcs
717945	cm 45 x cm 55	cm 45 x cm 49	10 pz/pcs
717950	cm 50 x cm 90	cm 50 x cm 84	10 pz/pcs
717960	cm 60 x cm 52	cm 60 x cm 46	10 pz/pcs
717970	cm 70 x cm 80	cm 70 x cm 74	10 pz/pcs
717980	cm 80 x cm 60	cm 80 x cm 54	10 pz/pcs
717990	cm 90 x cm 80	cm 90 x cm 74	10 pz/pcs

breathable



conformable



hypo-allergenic



latex free



# Euroderm incise iodine

## Outcomes efficacy verification test

stock	solutions	contact period (minutes)	stock UFC/ml	outcomes UFC/unit
staphylococcus aureus ATCC 25923	± 10.000 UFC/ml	15 30 60 to 300	11000	<10

- Reduction rate after 15 minutes: 100%

stock	solutions	contact period (minutes)	stock UFC/ml	outcomes UFC/unit
pseudomonas aeruginosa ATCC 27853	± 10.000 UFC/ml	15 30 60 to 300	10200	1000 80 <10

- Reduction rate after 15 minutes: 90,20%
- Reduction rate after 30 minutes: 99,20%T
- Reduction rate after 60 minutes: 100%

stock	solutions	contact period (minutes)	stock UFC/ml	outcomes UFC/unit
candida albicans ATCC 10231	± 10.000 UFC/ml	15 30 60 to 300	11000	2000 <10 <10

- Reduction rate after 15 minutes: 81,80%
- Reduction rate after 30 minutes: 100%

### Bactericide effect

